


MODEL NO.
SA-08

ENGINEERING

ENGINEERING

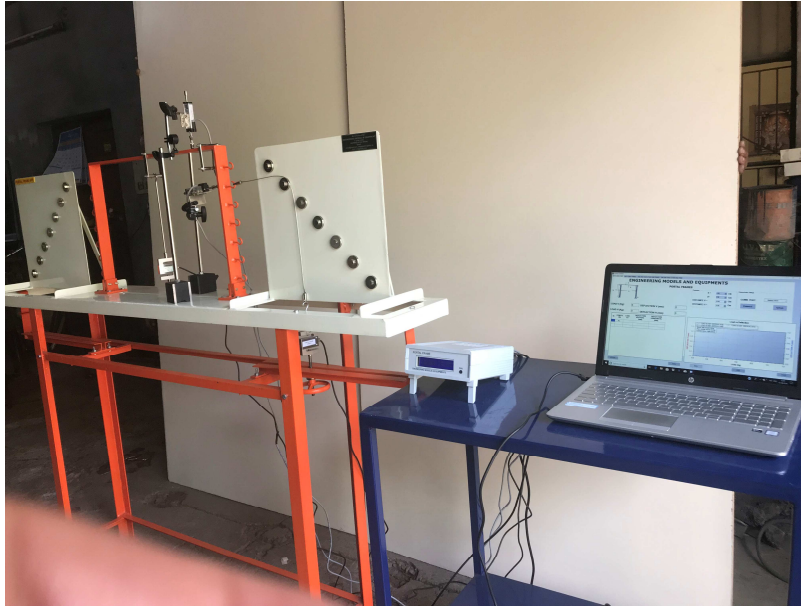

PORTAL FRAME APPARATUS (DIGITAL)

Description:

Portal frame is made up of M.S. plate of rectangular section of 8mm thick x 40cm wide. Frame is provided with a provision to achieve different end conditions viz. hinged, roller & fixed. The size of portal will be 40cm x 60cm. Portal is also having a provision for horizontal loading at different positions. Apparatus is supplied complete with a supporting stand and a set of weights.

Two LVDT with 25mm travel (with magnetic base) and two Load Cell of 10 Kg. capacities with a digital Indicator is supplied with the apparatus. Digital Indicator is capable to transfer the data through USB to the computer are supplied with the apparatus. Educational software also supplied with the set up to make calculation and complete results/data.

EDUCATIONAL SYSTEM YOU CAN BUY FROM:

ENGINEERING MODELS & EQUIPMENT

Regd. Office:

Mahavir Jain Market
7 - Civil Lines, P.B. No. - 13
Roorkee - 2476671 (U.K.)

Email: emodrke@gmail.com, enggmodels@bsnl.in

Factory & Works:

Adarsh Nagar
New Hardwar Road
Roorkee - 247667 (U.K.)
Phones: 01332-273721
Website: enggmod.com

AN ISO 9001:2000 & 14001:2004 COMPANY